

# HIRSCHMANN®



# ceramus®

HiClass in dispensing

# ceramus®

## Uncompromising reliability

Years of experience gained by Hirschmann in liquid handling are reflected in the perfection embodied in ceramus® and absolute attention to detail. The graduation gauge and digital volume configuration engagement guarantee consistent precision and volume dispensing and ensure reproducible results – absolutely independent of the user.

The dispensing piston made of high-purity aluminium oxide ceramic and the DURAN® glass\* cylinder guarantee a maximum level of chemical and thermal resistance.

This ensures continuously reliable dimensional stability of the parts which determine the volume and superior user comfort. The piston fits with absolute precision, can be operated easily with a minimum of effort and retains its functional reliability throughout the entire period of use.

The high-quality materials in all parts conveying medium enable the use of ceramus® within a temperature range from 4 °C to 50 °C, and the unit can be steam sterilised at 121 °C.

## Overview of details

- Piston made of high-purity aluminium oxide ceramic
- Fixed and variable volumes from 0.2 ml to 60 ml
- Complete unit can be sterilised at 121 °C
- Free rotation
- Media recirculating system for simple and rapid venting without reagent loss
- User-independent, reproducible volume dispensing
- Digital volume configuration engagement, calibrated for a lifetime
- Robust ejection unit
- Threaded suction hose fitting
- Optimum chemical and thermal resistance of parts conveying medium
- Quality certificate with declaration of conformity acc. to DIN 12600
- Discharge valves and valve seats consisting of high-purity aluminium oxide ceramic
- Valve springs from Hastelloy

\* DURAN® is a registered trademark of the DURAN Group GmbH, Wertheim.



**ceramus®**

## Absolute perfection down to the last detail

### consistent

Ceramic piston dimensional stability and functional reliability throughout the entire period of use

### high precision

Consistent volume dispensing and absolute reproducibility through graduation gauge and digital volume configuration engagement, calibrated for a lifetime

### ergonomic

A firm hold for control and simple, rapid, consistently exact volume configuration

### professional

User-independent  
consistent volume dispensing

### reliable

Quality certificate with declaration of conformity acc. to DIN 12600

### user-friendly

Valves can be unscrewed in seconds

### comprehensive

From 0.2 ml to 60 ml, fixed and variable

### Superior quality

Optimum chemical and thermal resistance of parts conveying media

### environmentally friendly

No reagent loss, thanks to user-friendly media recirculation system



Uncompromisingly reliable

# Precise and consistent

ceramus® is distinguished by a unique combination of a ceramic piston, graduation gauge and digital volume configuration engagement.

The chemical and thermal resistance of all parts conveying medium enables use of ceramus® for dispensing of practically all media within an extremely broad temperature range from 4 °C to 50°C. Functional reliability and precision remain absolutely guaranteed. The complete unit can be steam sterilised at 121 °C.



## Perfect ergonomics

Attention to detail: the ergonomically designed handle supports rapid, reliable and non-fatiguing working, both when heating the medium and dispensing.




## Precision instead of sense of proportion

Digital volume configuration engagement ensures immediate, consistent and absolutely reliable adjustment of the desired volume.



## Absolute reproducibility

The graduation gauge ensures absolutely reproducible work. A volume originally selected which is newly configured is always identical to the previous setting – with no need for subsequent calibration.

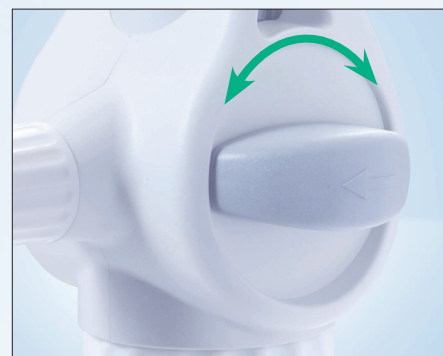

## User-friendly recirculation system

ceramus® is equipped with a unique comfortable media recirculation system. Venting is realised simply, reliably and rapidly with a single rotation – without reagent loss.




## Flexible and reliable

ceramus® can be screwed directly onto the bottle. The advantage: the labels or labelled side of the bottle remains visible to the user if, for example, the discharge is turned away from the user during dispensing of aggressive media. This means that you can keep an eye on everything.




### Precision work with suction hose

The robust suction hose is not pushed on, but screwed firmly to the unit. This reliably prevents the formation of air bubbles during suctioning.

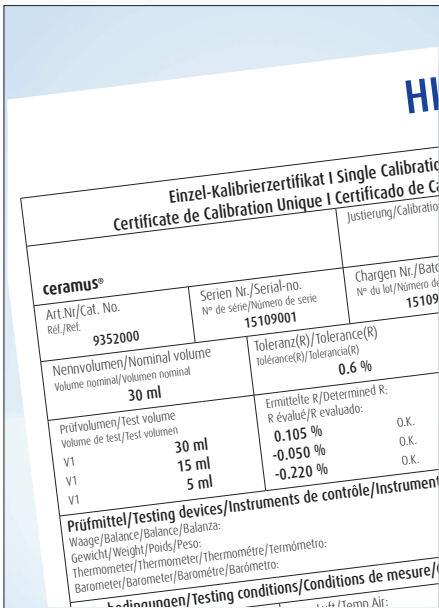

### Standard feature: quality certificate

Each ceramus® is individually inspected for precision. The quality certificate documents five tested volumes, the serial number and conformity pursuant to the German Calibration Regulation (Eichordnung) and can therefore be filed directly in quality documentation conforming to DIN EN ISO 9000. This ensures reliability.

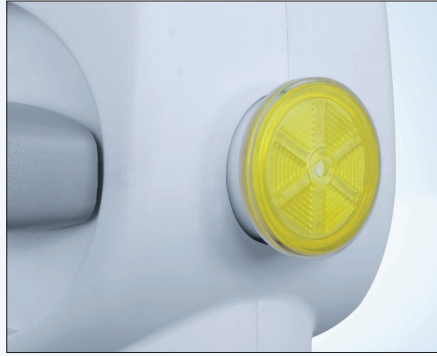

### Universal Luer connector

A universal Luer connector is fitted to the rear of the unit to enable connection of the dry tube, filters and much more.



### Ceramic piston made of extremely hard, high-purity aluminium oxide ceramic

The ceramic piston is made of extremely hard, high-purity aluminium oxide ceramic and is fitted into the DURAN® glass cylinder with micrometre precision for an absolutely tight seal. This achieves impressive ease of guidance and a maximum degree of functional reliability. Abrasion, serration, gushing or the detaching of a layer are prevented.

# ceramus® classic

## The classic unit for rapid dispensing with reliable results

ceramus® classic is the proven basic unit for manual dispensing. Superior-quality materials, simple handling and well-engineered user comfort make laboratory routines easier.

The graduation gauge and digital volume configuration engagement also guarantee consistent precision and volume dispensing and ensure reproducible results in the case of the ceramus® classic – absolutely independent of the user. The dispensing piston made of high-purity aluminium oxide ceramic guarantees a maximum level of chemical and thermal resistance. The complete unit can be steam sterilised at 121 °C.

The ceramus® classic FIX is also available in six different volumes for dispensing with fixed volume configurations.

### Overview of details

- Piston made of high-purity aluminium oxide ceramic
- For fixed and variable volumes from 0.2 ml to 60 ml
- Complete unit can be steam sterilised at 121 °C
- Threaded suction hose fitting
- User-independent, reproducible volume dispensing
- Digital volume configuration engagement, calibrated for a lifetime
- Optimum chemical and thermal resistance of parts conveying medium
- Individual quality certificate with declaration of conformity acc. to DIN 12600
- Discharge valves and valve seats consisting of high-purity aluminium oxide ceramic
- Valve springs from Hastelloy



## **ceramus® HF**

**Specially developed  
for hydrofluoric acid or  
other aggressive media**

A special version designed for dispensing very aggressive media, the ceramus® HF, is also available. The piston, cylinder, valve seats and balls are made of high-purity aluminium oxide ceramic, while the valve springs are made of platinum-iridium.

The volume can be variably adjusted from 2 ml to 10 ml.

## **EM-dispenser PP**

The EM Dispenser PP with parts for conveying media made of polypropylene supplements the ceramus® family as a light version. With a precision ground glass piston and cylinder made of DURAN® glass, it is ideal for non-problematic watery and organic media which do not attack polypropylene.

ceramus® family



| Order number | Article designation   | ceramus® | ceramus® classic up to 10 ml and EM Dispenser PP 2-10 ml | ceramus® classic from 25 ml and EM Dispenser PP 10-60 ml | ceramus® HF |
|--------------|-----------------------|----------|----------------------------------------------------------|----------------------------------------------------------|-------------|
|              | Threaded adapters     |          |                                                          |                                                          |             |
| 931 35 23    | PP A 45 to A 28       | ✓        |                                                          |                                                          | ✓           |
| 931 35 24    | PP A 45 to A 38/430   | ✓        |                                                          |                                                          | ✓           |
| 931 36 21    | ETFE A 45 to A 32     | ✓        |                                                          |                                                          | ✓           |
| 931 36 23    | ETFE A 45 to A 28     | ✓        |                                                          |                                                          | ✓           |
| 931 36 24    | ETFE A 45 to A 38/430 | ✓        |                                                          |                                                          | ✓           |
| 931 36 25    | ETFE A 45 to A 38     | ✓        |                                                          |                                                          | ✓           |
| 931 36 27    | ETFE A 45 to S 40     | ✓        |                                                          |                                                          | ✓           |
| 931 35 02    | PP A 32 to A 25       |          | ✓                                                        | ✓                                                        |             |
| 931 35 05    | PP A 32 to A 38       |          | ✓                                                        | ✓                                                        |             |
| 931 36 06    | ETFE A 32 to A 45     |          | ✓                                                        | ✓                                                        |             |
| 931 36 07    | ETFE A 32 to S 40     |          | ✓                                                        | ✓                                                        |             |

|           |                                                                                                                                           |   |   |   |   |
|-----------|-------------------------------------------------------------------------------------------------------------------------------------------|---|---|---|---|
|           | Valves                                                                                                                                    |   |   |   |   |
| 931 60 28 | Discharge valve, valve seat/ball made of high purity aluminium oxide ceramic ( $Al_2O_3$ 99,7%), spring made of platinum-iridium (HF)     | ✓ |   |   | ✓ |
| 931 60 10 | Discharge valve, valve seat/ball made of high purity aluminium oxide ceramic ( $Al_2O_3$ 99,7%), spring made of Hastelloy                 |   | ✓ |   |   |
| 931 60 15 | Discharge valve, valve seat/ball made of high purity aluminium oxide ceramic ( $Al_2O_3$ 99,7%), spring made of Hastelloy                 | ✓ |   | ✓ |   |
| 931 65 28 | Suction valve, valve seat/ball made of high purity aluminium oxide ceramic ( $Al_2O_3$ 99,7%), spring made of platinum-iridium (HF)       | ✓ |   |   | ✓ |
| 931 65 10 | Suction valve, valve seat/ball made of high purity aluminium oxide ceramic ( $Al_2O_3$ 99,7%), spring made of Hastelloy                   |   | ✓ |   |   |
| 931 65 15 | Suction valve, valve seat/ball made of high purity aluminium oxide ceramic ( $Al_2O_3$ 99,7%), spring made of Hastelloy                   | ✓ |   | ✓ |   |
| 931 66 26 | Recirculation valve, valve seat/ball made of high purity aluminium oxide ceramic ( $Al_2O_3$ 99,7%), spring made of platinum-iridium (HF) | ✓ |   |   | ✓ |
| 931 66 15 | Recirculation valve, valve seat/ball made of high purity aluminium oxide ceramic ( $Al_2O_3$ 99,7%), spring made of Hastelloy             | ✓ |   |   |   |

|           |                                                                        |   |   |   |   |
|-----------|------------------------------------------------------------------------|---|---|---|---|
|           | Ejection units                                                         |   |   |   |   |
| 931 30 15 | Titration ejection unit                                                | ✓ |   |   | ✓ |
| 956 30 04 | Ejection unit, titration, variable, with hose clip, 400 mm, FEP        | ✓ |   |   | ✓ |
| 931 30 26 | Ejection unit, flexible, spiral shaped                                 | ✓ |   |   | ✓ |
| 931 30 27 | Ejection unit, flexible, spiral shaped                                 | ✓ |   |   | ✓ |
| 931 30 28 | Ejection unit, with Luer lock connector                                | ✓ |   |   | ✓ |
| 931 30 29 | Ejection unit, with Luer lock connector                                | ✓ |   |   | ✓ |
| 931 30 07 | Ejection unit FEP, flexible, spiral shaped, working length approx. 1 m |   |   | ✓ |   |
| 931 30 08 | Ejection unit, with Luer lock connector                                |   |   | ✓ |   |
| 931 30 09 | Ejection unit, with Luer lock connector                                |   |   | ✓ |   |
| 931 30 09 | Ejection unit FEP, flexible, spiral shaped, working length approx. 1 m |   | ✓ |   |   |

|           |                                 |   |   |   |   |
|-----------|---------------------------------|---|---|---|---|
|           | Suction/ Recirculation hoses    |   |   |   |   |
| 931 50 07 | Suction hose with screw fitting | ✓ |   | ✓ | ✓ |
| 931 50 07 | Suction hose FEP with union nut |   | ✓ |   |   |

|           |                                                     |   |  |  |   |
|-----------|-----------------------------------------------------|---|--|--|---|
|           | Miscellaneous                                       |   |  |  |   |
| 990 37 04 | Round filter set (with plug connection on one side) | ✓ |  |  | ✓ |

## Newly developed accessory range for maximum adaptability



### Base support instead of bottle attachment

The base support enables suctioning of the medium directly from a separate storage container using a hose.



### Stand clip for fixing the basic unit

The stand clip also facilitates use of ceramus® in a stand assembly.



### Adapter with NS 12/21 filling tube for funnel

The adapter with standard ground 12/21 filling tube enables problem-free filling of medium through a funnel - without unscrewing the basic unit.



### Ground adapter for direct reagent extraction from a ground bottle

The NS 29/32 ground adapter enables secure use of the basic unit on standard ground bottles.

Accessory range



HIRSCHMANN®

10-60 ml

div. in 0,5 ml Ex 20°C

DE-M16

ceramus®

HIRSCHMANN

5 ml

div. in 0,1 ml

ceramus

200

ceramus® HiClass in dispensing

| Performance data    |                                                                                                                                                                                                                                                                                        |
|---------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Volume selection    | With micro-engagement adjusting ring, precise and reproducible, thanks to mechanical digital principle (with fixed preadjustment in defined graduations); individual graduations with 0.05 to 1 ml resolution (depending on model)                                                     |
| Volume accuracy     | Precise and reproducible with absolute accuracy (R) of 0.6% and standard deviation (CV) of 0.2%                                                                                                                                                                                        |
| Piston              | Mainly inert and wear-free ceramic consisting of 99.7% pure Al <sub>2</sub> O <sub>3</sub> (EM Dispenser PP made of Duran®)                                                                                                                                                            |
| Cylinder            | Mainly inert and wear-free design made of calibrated DURAN® glass (HF: ceramic)                                                                                                                                                                                                        |
| Valves              | Mainly inert and wear-free design, ball and seat made of ceramic (99.7% pure Al <sub>2</sub> O <sub>3</sub> ), ejection valve springs made of Hastelloy (HF: platinum-iridium). The valves can be easily removed and replaced using the fitting wrench (included in scope of delivery) |
| Plastic parts       | All parts which come into contact with the medium are made of highly resistant fluoroplastics: Valve housing made of ECTFE, suction and ejection hose made of FEP (EM Dispenser PP: valve block made of PP)                                                                            |
| Sterilisability     | Can be steam sterilised at 121°C, complete including ejection unit                                                                                                                                                                                                                     |
| Breakage resistance | Completely enclosed protective plastic jacket with clear inspection window                                                                                                                                                                                                             |
| Test protocol       | Each instrument is individually tested, certified and bears a unique instrument number                                                                                                                                                                                                 |


| Scope of delivery                                                                                                                                                                                                                                                                                                                                     |
|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <ul style="list-style-type: none"><li>ceramus®, suitable for A 45 thread; ceramus® classic, suitable for A 32 thread</li><li>Suction hose, fitting wrench, individual quality certificate, ceramus® additionally with recirculation hose</li><li>Threaded adapters A 32, A 38 and S 40, ceramus® classic threaded adapters A 28, S 40, A 45</li></ul> |

| Ordering data               |             |                 |                |                       |                                          |
|-----------------------------|-------------|-----------------|----------------|-----------------------|------------------------------------------|
|                             | Volume (ml) | Graduation (ml) | Precision (R%) | Reproducibility (CV%) | Order no.<br>ceramus®   ceramus® classic |
| ceramus®   ceramus® classic | 0.2 - 1.0   | 0.05            | 0.6            | 0.2                   | 931 20 00   931 00 00                    |
| ceramus®   ceramus® classic | 0.4 - 2.0   | 0.05            | 0.6            | 0.2                   | 932 20 00   932 00 00                    |
| ceramus®   ceramus® classic | 1.0 - 5.0   | 0.10            | 0.6            | 0.2                   | 933 20 00   933 00 00                    |
| ceramus®   ceramus® classic | 2.0 - 10.0  | 0.25            | 0.6            | 0.2                   | 934 20 00   934 00 00                    |
| ceramus®   ceramus® classic | 5.0 - 30.0  | 0.50            | 0.6            | 0.2                   | 935 20 00   935 00 00                    |
| ceramus®   ceramus® classic | 10.0 - 60.0 | 1.00            | 0.6            | 0.2                   | 936 20 00   936 00 00                    |
| ceramus®   ceramus® classic | 1 FIX       | -               | 0.6            | 0.2                   | 931 21 00   931 11 00                    |
| ceramus®   ceramus® classic | 2 FIX       | -               | 0.6            | 0.2                   | 932 21 00   932 11 00                    |
| ceramus®   ceramus® classic | 5 FIX       | -               | 0.6            | 0.2                   | 933 21 00   933 11 00                    |
| ceramus®   ceramus® classic | 10 FIX      | -               | 0.6            | 0.2                   | 934 21 00   934 11 00                    |
| ceramus®   ceramus® classic | 25 FIX      | -               | 0.6            | 0.2                   | 935 21 00   935 11 00                    |
| ceramus®   ceramus® classic | 50 FIX      | -               | 0.6            | 0.2                   | 936 21 00   936 11 00                    |
| ceramus® HF                 | 2.0 - 10.0  | 0.25            | 0.6            | 0.2                   | 934 23 00                                |
| EM-dispenser PP             | 2.0 - 10.0  | 0.25            | 0.6            | 0.2                   | 934 02 00                                |
| EM-dispenser PP             | 10.0 - 60.0 | 1.00            | 0.6            | 0.2                   | 936 02 00                                |

ceramus® HiClass dispensing

# HIRSCHMANN®

Hirschmann Laborgeräte GmbH & Co. KG  
Hauptstraße 7-15 • 74246 Eberstadt Germany  
Fon +49 7134 511-0 • Fax +49 7134 511-990  
[www.hirschmannlab.com](http://www.hirschmannlab.com) • [info@hirschmannlab.com](mailto:info@hirschmannlab.com)

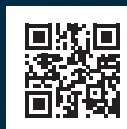

X2100450100