


BIOHIT

Innovating for Health

eLINE® Dispenser

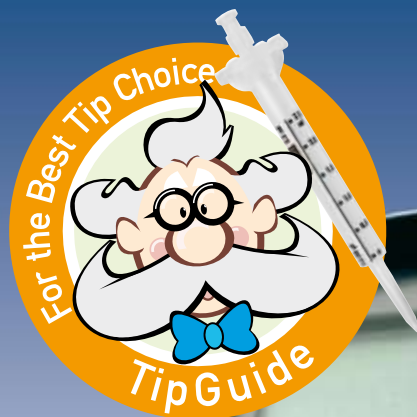

A single innovation makes all the difference

eLINE electronic dispenser hands you a familiar but unique approach to a new, more comfortable world of multiple dispensing.

Its unique TipGuide selects and displays the most suitable tip size for you - saving time and eliminating the possibility of an incorrect tip selection. The automated liquid handling operations offer accurate and repetitive dispensing also with viscous, infectious and high vapour pressure liquids. Additionally, a smart feature included in the whole eLINE product range is the unique electronic tip ejection, which reduces the risk of RSI (Repetitive Strain Injury).

eLINE Dispenser is your most convenient and ergonomic choice.





eLINE Dispenser offers

- TipGuide – automatic guide for correct tip selection
- Outstanding ergonomics – dispensing with minimum force
- Unique electronic feather light tip ejection
- Comprehensive range of liquid handling protocols with easy programming
- Contamination-free dispensing with positive displacement principle
- Volume range from 1 µl to 50 ml
- 5 speeds for aspiration and dispensing
- Direct charging option enabling dispensing while charging
- Compatibility with eLINE charging stand and carousel
- Compatibility with other dispenser tip brands

eLINE Dispenser operating modes

- Multiple Dispensing (d)
- Pipetting (P)
- Diluting (dd)
- Sequential Dispensing (Sd)
- Automatic Multiple Dispensing (Ad)
- Multi-Aspirating (SA)
- Custom mode (CST)

eLINE Dispenser Ordering details

Cat. No.	Item
73070X	Biohit eLINE Dispenser
73098X	Biohit eLINE Charging Stand, with AC adaptor
73099X	eLINE Charging Carousel for 4 pipettors, with AC adaptor
731001	eLINE replacement battery

X: 0 =Without AC-adaptor, 1 = With AC-adaptor, Euro; 2 = With AC-adaptor, U.S.; 3 = With AC-adaptor, U.K.; 4 = With AC-adaptor, Jpn; 5 = With AC-adaptor, Aus
Biohit has several worldwide patents and patents pending for its products.

Distributed by:



Save your time, let the TipGuide make the tip choice for you

- No time consuming volume calculations
- Just set the desired volume and number of dispensings, the TipGuide selects and displays the suitable tip size

Clear LCD-display with charging indicator

Easy-to-use keypad enables fast programming

Unique electronic one-touch tip ejection

Electronic start button for all pipetting operations

Easy & light tip insertion

Direct charging option with AC-adaptor with immediate dispensing possibility

Both left- and right-handed operation



Biohit Dispenser Tips Ordering details

Cat. No.	Volume [ml]	Type	Qty/Unit
792017	0.1	Non sterile	100
792018	0.2	Non sterile	100
792019	0.5	Non sterile	100
792020	1.0	Non sterile	100
792021	2.5	Non sterile	100
792022	5.0	Non sterile	100
792023	10.0	Non sterile	100
792024	25.0*	Non sterile	25*
792025	50.0*	Non sterile	25*
792036	Adapter for 25 and 50 ml tips, 3 pcs./unit		

* includes 1 adapter/unit



For more information and a list of all Biohit distributors visit www.biohit.com